



Using Oxytocin in the Uniject™* Injection System (10 IU in 1 ml)


- 1** Check the time-temperature indicator and confirm the oxytocin is ok to use. If not, discard and get a new Uniject containing oxytocin.

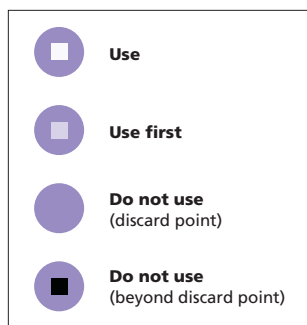


- 2** Open the foil pouch and remove the Uniject.

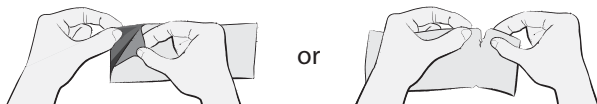



- 3** Activate the Uniject.
- Hold the Uniject by the port with the forefinger and thumb.
 - With a firm, rapid movement push the needle shield into the port.
 - Continue to push firmly until the gap between the needle shield and port is closed.

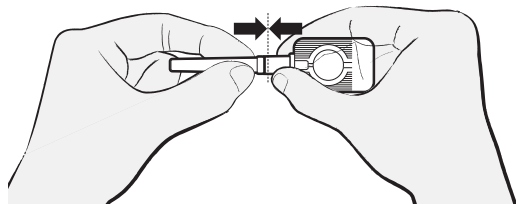


- 4** Remove the needle cap.

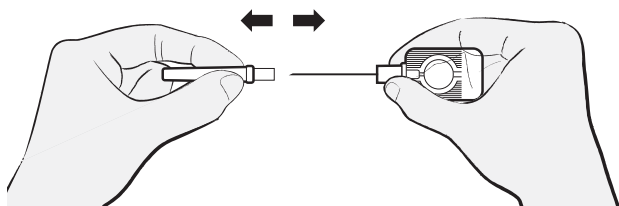


- 5** Continue to hold the Uniject by the port and insert the needle into the patient.

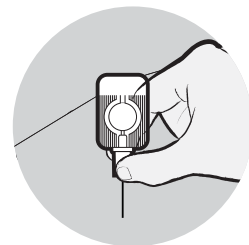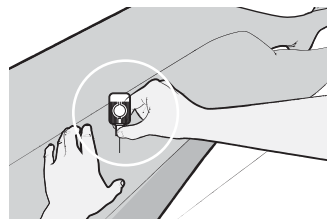


- 6** Squeeze the reservoir firmly to inject the oxytocin. After the reservoir completely collapses, remove the Uniject.

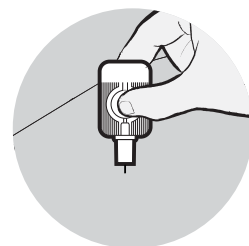

- 7** Do **not** re-cap the used Uniject.

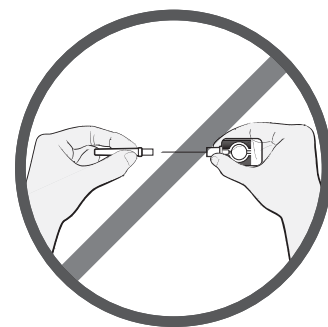

- 8** Discard the used Uniject according to established medical waste disposal procedures.

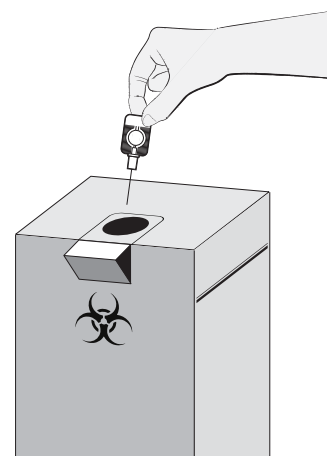

* Uniject is a trademark of BD.

Copyright © 2011, Program for Appropriate Technology in Health (PATH). The material in this document may be freely used for educational or noncommercial purposes, provided that the material is accompanied by an acknowledgment line.