


Reading time-temperature indicators with the Uniject™ injection system


Worksheet

Complete the spaces below with one of the following: **USE**, **USE FIRST**, or **DON'T USE**.


Inner square is lighter than outer circle.

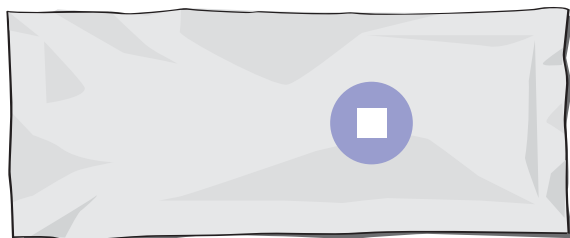


Inner square is lighter than outer circle.

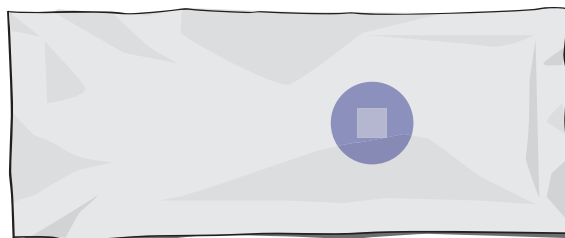

Inner square matches outer circle.

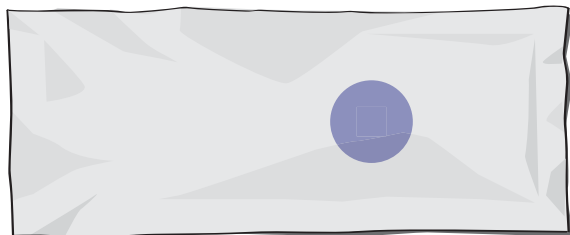

Inner square is darker than outer circle.

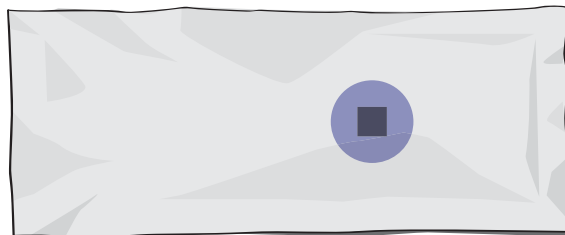


Copyright © 2011, Program for Appropriate Technology in Health (PATH). The material in this document may be freely used for educational or noncommercial purposes, provided that the material is accompanied by an acknowledgment line.

Reading time-temperature indicators with the Uniject™ injection system

Answers


USE

Inner square is lighter than outer circle.

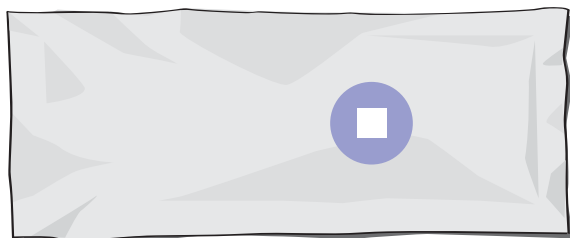

USE FIRST


Inner square is lighter than outer circle.

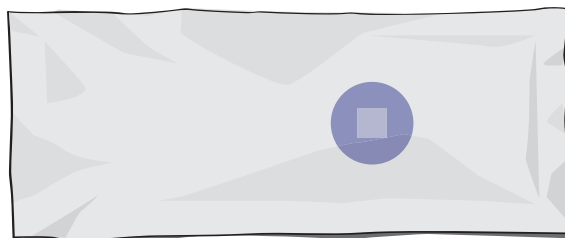

DON'T USE

Discard point!


Inner square matches outer circle.

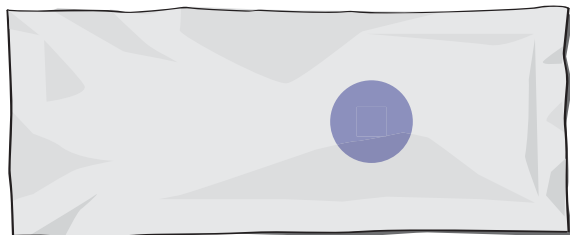

DON'T USE

Beyond discard point!

Inner square is darker than outer circle.

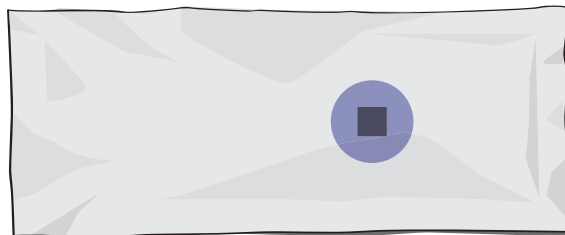

Copyright © 2011, Program for Appropriate Technology in Health (PATH). The material in this document may be freely used for educational or noncommercial purposes, provided that the material is accompanied by an acknowledgment line.